

Inertgas-Systems

Gloveport oval-long Operating Manual



Original instructions

Legal information

Operating Manual

Gloveport oval-long

M. Braun Inertgas-Systeme GmbH
Dieselstr. 31
85748 Garching
Germany

Website: www.mbraun.de
Phone: +49 (0)89 32669-0
Fax: +49 (0)89 32669-235
E-Mail: service@mbraun.de

1 Maintenance

1.1 Replacing a glove




This manual only describes the glove replacement on oval-long gloveports. For all other information regarding gloves (maintenance intervals, user information etc.) refer to the main system documentation.

Mounting the gloveport cover (optional)

To prevent contaminating the glovebox atmosphere by ambient air during a glove replacement, the gloveport can be sealed from inside the box with an optionally available gloveport cover.

Prerequisite

- > The O-rings of the gloveport cover are undamaged and free of dust.
- > The gloveport cover is inside the glovebox.
- > The glovebox is set to slight overpressure (approx. 1 ... 2 mbar).

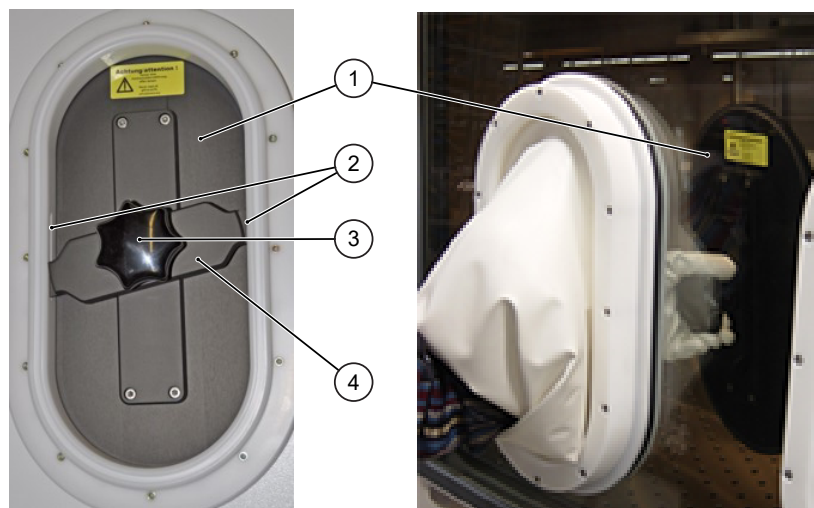

- ▶ Turn the star knob (3) counterclockwise as far as possible.
- ▶ Grab the gloveport cover (1) with the glove by the star knob (3) and pull it outward to attach it to the gloveport.
- ▶ Make sure that the gloveport cover (1) is properly seated in the gloveport.
- ▶ Begin to turn the star knob (3) clockwise.
- ⇒ The lock bar (4) should now snap into the notches (2).
- ▶ If the lock bar (4) did not snap into the notches (2) completely, turn it into the notches (2) by hand.
- ▶ Turn the star knob (3) clockwise as far as possible to fasten the gloveport cover (1) tightly.

Removing the glove

Required materials


- Fitting spare glove (type, size etc.)
- Screwdriver (hexagon socket, M5, torque approx. 120 N/cm)
- Disposal container for the glove (if necessary)

Prerequisites

- > The gloveport is tightly sealed with the gloveport cover (if applicable).
- > If no gloveport cover is used, the glovebox must be set to slight overpressure (approx. 1 ... 2 mbar).



- Remove the screws of the clamp ring.

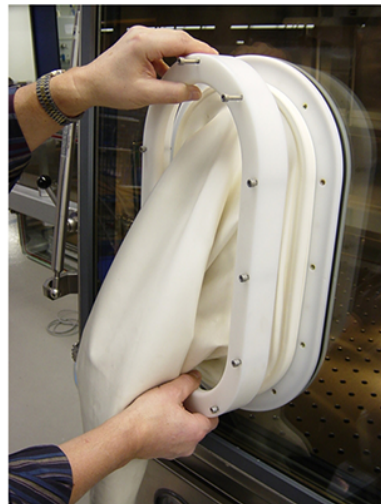

- Remove the clamp ring.




- Remove the used glove and put it into the disposal container (if necessary).

Mounting the new glove



- Roll up the new glove to remove as much air as possible.

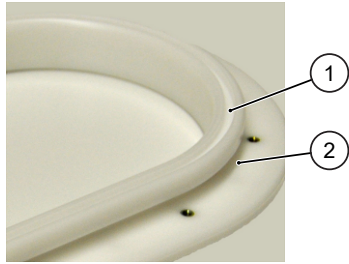

- ▶ Check the O-ring in the O-ring notch (1) and replace it if necessary.
- ▶ Slip the shaft of the new glove over the glove notch (2).



- ▶ Make sure that the shaft fits smooth and is without wrinkles.



- ▶ Place the rolled glove in the gloveport opening.



- ▶ Attach the new clamping ring, avoiding pinching or damaging the new glove.



- ▶ Put the screws back in place and hand-tighten them.
- ▶ Tighten all screws with a tightening torque of 2 Nm.
- ▶ Afterwards tighten all screws with a final torque of 3 Nm.
- ▶ Turn the star knob of the gloveport cover counterclockwise and remove it slowly (if applicable).
- ▶ Set the glovebox to normal working pressure.
- ▶ Purge the glovebox if necessary (refer to the main system documentation).

2 Spare parts

Gloves




No.	Item	Description	Item no.
–	CSM glove, oval, size 9.5 / L, 2 pieces, ambidextrous	D = 300, 800 x 0,4 mm	2127050



No.	Item	Description	Item no.
–	Antistatic butyl glove, oval, size 9.5 / L, 2 pieces, ambidextrous	D = 300, 800 x 0,4 mm	3005008
–	Antistatic butyl glove, oval, size 9.5 / L, 2 pieces, ambidextrous	D = 250, 800 x 0,4 mm	2608918

Gloveport

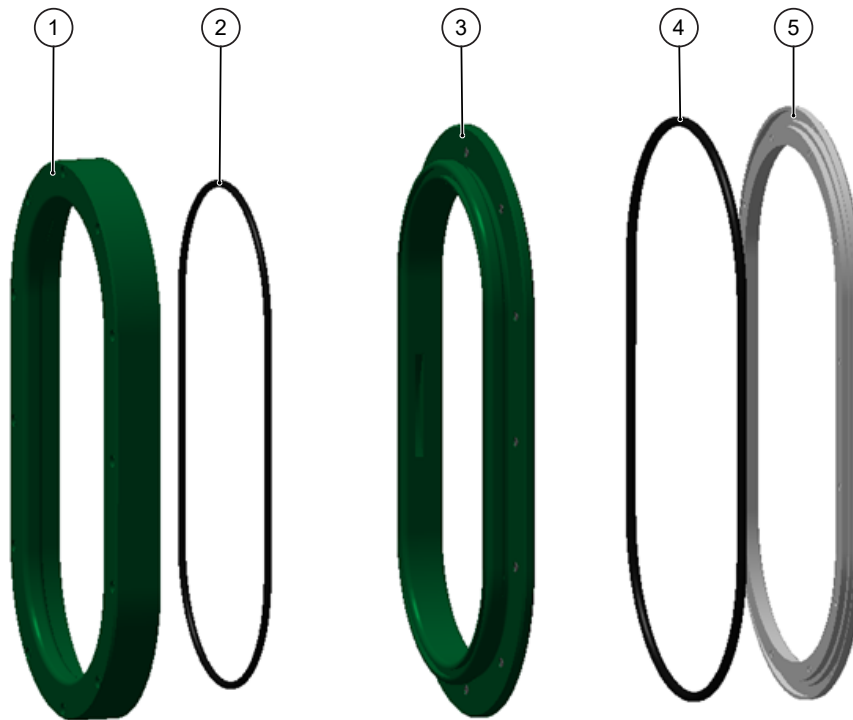

No.	Item	Description	Item no.
–	Gloveport complete, POM-C	with O-Rings	1515405
	Gloveport complete, POM-C-ELS, antistatic		7098950
1	Clamping ring, POM-C	444 x 244 x 38	7041046
	Clamping ring, POM-C-ELS, antistatic		7098952
2	O-rings for gloveport, outside	320 x 6	2400382r
3	Outer ring, POM-C	470 x 270x 34	7041045
	Outer ring, POM-C-ELS, antistatic		7098951
4	O-rings for gloveport, inside	374 x 8	2400383r
5	Inner ring, POM-C	470 x 270 x 13	7041047
	Inner ring, POM-C-ELS, antistatic		7098953

Gloveport cover



No.	Item	Description	Item no.
-	MB-HVD-Internal gloveport cover	D = 300	1515203

